



Model: AIRPAP ASV Serial: SV47259999 SW: 1.06

06 May 2026 – 06 May 2026 | 3 sessions | 3h 05m

Patient Information Name: **Gökhan Gültekin** Gender: **Male** Age: **46 years** Height: **172 cm** Weight: **74,0 kg**

MULTI-RECORD SUMMARY

06 May 2026 – 06 May 2026

3 sessions

SLEEP QUALITY SCORE

EXCELLENT THERAPY

A

RECORD START

06 May 2026 10:04:12

RECORD END

06 May 2026 16:09:08

TOTAL DURATION

3h 05m

AHI

Apnea-Hypopnea

Ø 0,8

↓0,0 ↑2,3

• Normal (<5)

OAI

Obstructive Apn

Ø 0,6

↓0,0 ↑1,7

• Normal (<5)

CAI

Central Apnea

Ø 0,0

↓0,0 ↑0,0

• Normal (<5)

HI

Hypopnea

Ø 0,2

↓0,0 ↑0,6

• Normal (<5)

SI


Snore

Ø 0,0

↓0,0 ↑0,0

• Normal (<5)

AHI per Session

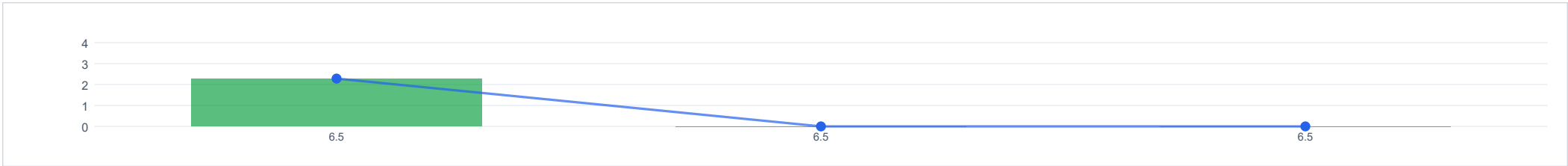

Session Overview

#	Session Start	Session Stop	Duration	AHI	OAI	HI	Closed
1	06 May 2026 10:04	06 May 2026 11:49	1h 45m	2,3	1,7	0,6	4
2	06 May 2026 14:48	06 May 2026 15:26	0h 38m	0,0	0,0	0,0	0
3	06 May 2026 15:27	06 May 2026 16:09	0h 42m	0,0	0,0	0,0	0



Start Parameters by Mode — values across sessions					
Mode	Idx	Parameter Name	Min	Max	↑
BiLEVEL-T	17	Inspiration Pressure	15,0	15,0	
	18	Expiration Pressure	5,0	5,0	
	19	Inspiration Slope	1,0	1,0	
	20	Expiration Slope	4,0	4,0	
	21	Frequency	12,0	12,0	
	22	Inspiration Duration	1,0	1,0	
	23	Assist Mode	0,0	0,0	
	24	Sigh Function	0,0	0,0	



Model: AIRPAP ASV Serial: SV47259999 SW: 1.06

06 May 2026 – 06 May 2026 | 3 sessions | 3h 05m

Page 3 / 4

Aggregate Statistics — All Sessions

Statistic	Min	Mean (Ø)	Max	Unit
AHI / Respiratory Indices				
AHI	0,0	0,8	2,3	events/h
OAI	0,0	0,6	1,7	events/h
CAI	0,0	0,0	0,0	events/h
HI	0,0	0,2	0,6	events/h
SI	0,0	0,0	0,0	events/h
Session Duration				
Session Duration	38,7	62,0	105,1	min
Event Counts				
Total Closed Events	0,0	1,3	4,0	events
APNEA	0,0	1,0	3,0	events
HIPOPNEA	0,0	0,3	1,0	events
CENTRAL APNEA	0,0	0,0	0,0	events
SNORE	0,0	0,0	0,0	events
Parameter Changes				
Param Changes / Session	0,0	2,3	5,0	changes
Mode Changes / Session	0,0	0,7	2,0	changes



Model: AIRPAP ASV Serial: SV47259999 SW: 1.06

06 May 2026 – 06 May 2026 | 3 sessions | 3h 05m

Page 4 / 4

Sleep Events — All Sessions

Session	Name	Start	End	Duration	Closed
#1	APNEA	10:28:12-6.5.2026	10:28:13-6.5.2026	145,8s	Yes
#1	APNEA	10:38:51-6.5.2026	10:50:43-6.5.2026	73,3s	Yes
#1	HIPOPNEA	10:50:58-6.5.2026	10:50:58-6.5.2026	1,1s	Yes
#1	APNEA	10:50:59-6.5.2026	10:50:59-6.5.2026	1,1s	Yes

Param Changes — All Sessions

Session	DateTime	Mode	Changed Parameter	New Value
#1	10:9:57-6.5.2026	BiLEVEL-T	[17] Inspiration Pressure	20,0
#1	10:18:31-6.5.2026	BiLEVEL-T	[21] Frequency	20,0
#1	10:28:12-6.5.2026	BiLEVEL S	MODE CHANGE	
#1	10:30:40-6.5.2026	BiLEVEL S	[13] Inspiration Sensitivity	3,0
#1	11:3:20-6.5.2026	BiLEVEL-T	MODE CHANGE	
#1	11:17:45-6.5.2026	BiLEVEL-T	[17] Inspiration Pressure	15,0
#1	11:29:31-6.5.2026	BiLEVEL-T	[21] Frequency	12,0
#2	14:58:31-6.5.2026	BiLEVEL-T	[21] Frequency	20,0
#2	15:8:15-6.5.2026	BiLEVEL-T	[17] Inspiration Pressure	20,0